

## Table of contents

1. Overview
2. Contact
3. Disclaimer
4. Installation notes
5. Running the software
6. Basic operation



## 1. Overview

T2mapp is a software utility and graphical user interface to reconstruct T2 maps from multi-echo spin-echo and T2\* maps from multi-gradient-echo MRI acquired with an MR Solutions preclinical MRI system.

## 2. Contact

T2mapp and this manual are written by:

Gustav Strijkers  
Amsterdam UMC  
Biomedical Engineering and Physics  
Amsterdam, the Netherlands  
[g.j.striekers@amsterdamumc.nl](mailto:g.j.striekers@amsterdamumc.nl)

## 3. Disclaimer

The software has been extensively tested using mouse and rat data acquired with an MR Solutions preclinical 7.0T/24cm system.

However, this does not warrant the functions contained in the program will meet your requirements or that the operation of the program will be uninterrupted or error-free.

In case of questions or issues, please contact Gustav Strijkers or Bram Coolen.



### 4. Installation notes

#### Software download

Matlab source code and a Windows standalone version (using the free Matlab runtime engine) can be downloaded from github:

<https://github.com/Moby1971?tab=repositories>

#### Installation of the Windows standalone version

`MyAppInstaller_web.exe`

Will install the Matlab runtime engine and the T2mapp program.

#### Bart toolbox download

The advanced reconstruction options in P2ROUD require the Bart reconstruction toolbox, which can be downloaded from:

<https://mrirecon.github.io/bart/>



### Bart toolbox installation in OSX with MacPorts

- (1) Install Xcode from the Mac App Store

```
$ xcode-select --install
```

- (2) Install MacPorts (<http://www.macports.org/>)

It is recommended to install a newer version of gcc from MacPorts

- (3) Install packages

```
$ sudo port install fftw-3-single
$ sudo port install gcc10
$ sudo port install libpng
$ sudo port install openblas
```

- (4) Make

Replace 'gcc-mp-6' with 'gcc-mp-10' in 'makefile'

```
$ make
```

### Bart toolbox installation in OSX with HomeBrew (for Mac with Apple Silicon)

- (1) Install Xcode from the Mac App Store

```
$ xcode-select --install
```

- (2) Install HomeBrew

```
$ /bin/bash -c "$(curl -fsSL https://raw.githubusercontent.com/Homebrew/install/HEAD/install.sh)"
```

- (3) Install packages

```
$ brew install --cask gcc-arm-embedded
$ brew install libpng
$ brew install fftw
$ brew install openblas
```

- (3) Make

- Use the provided makefile for Apple Silicon

```
$ make
```



### Bart toolbox installation in Windows 10

#### (1) Install the Windows subsystem for Linux

Start a windows powershell as administrator and run the following command:

```
Enable-WindowsOptionalFeature -Online -FeatureName Microsoft-Windows-Subsystem-Linux
```

A system restart will be needed.

For more information, see:

<https://docs.microsoft.com/en-us/windows/wsl/install-win10>

#### (2) Download Ubuntu Linux

Open PowerShell

To see a list of available Linux distributions for download, enter:

```
wsl -l -o
```

Install a Linux distribution:

```
wsl --install -d <distribution name>
```

You will be asked to create a user account the first time you start a Linux command prompt.

#### (3) Upgrade the Linux distribution

```
$ sudo apt-get update
```

```
$ sudo apt-get dist-upgrade
```

#### (4) Install Bart prerequisites

```
$ sudo apt-get install make gcc gcc-10 libfftw3-dev liblapacke-dev  
libpng-dev libopenblas-dev gfortran
```

#### (5) Download Bart

```
$ wget https://github.com/mrirecon/bart/archive/v0.7.00.tar.gz
```

For WSL1 version 0.4.02 is recommended. Newer versions seem to compile but some functions which are required produce errors.

Recommended: Upgrade to WSL2 (see e.g. <https://pureinfotech.com/install-windows-subsystem-linux-2-windows-10/> on how to upgrade).

#### (6) Build Bart

```
$ tar xvfz v0.7.00.tar.gz
```

```
$ cd bart-0.7.00
```

```
$ make
```

```
$ sudo make PREFIX=/usr/local install
```



## 5. Running the software

### Running in Matlab 2021b

The T2mapp software can be started from its root directory from the command line.

```
>> T2mapp
```

Notes:

Additional licenses may be required.

```
>> license 'inuse'
distrib_computing_toolbox
image_toolbox
matlab
phased_array_system_toolbox
signal_blocks
signal_toolbox
```


### Running the Windows standalone

The Windows standalone version can be run from the start menu or the desktop icon.



## 6. Basic operation

The T2mapp program operates from a single window with 5 panels.

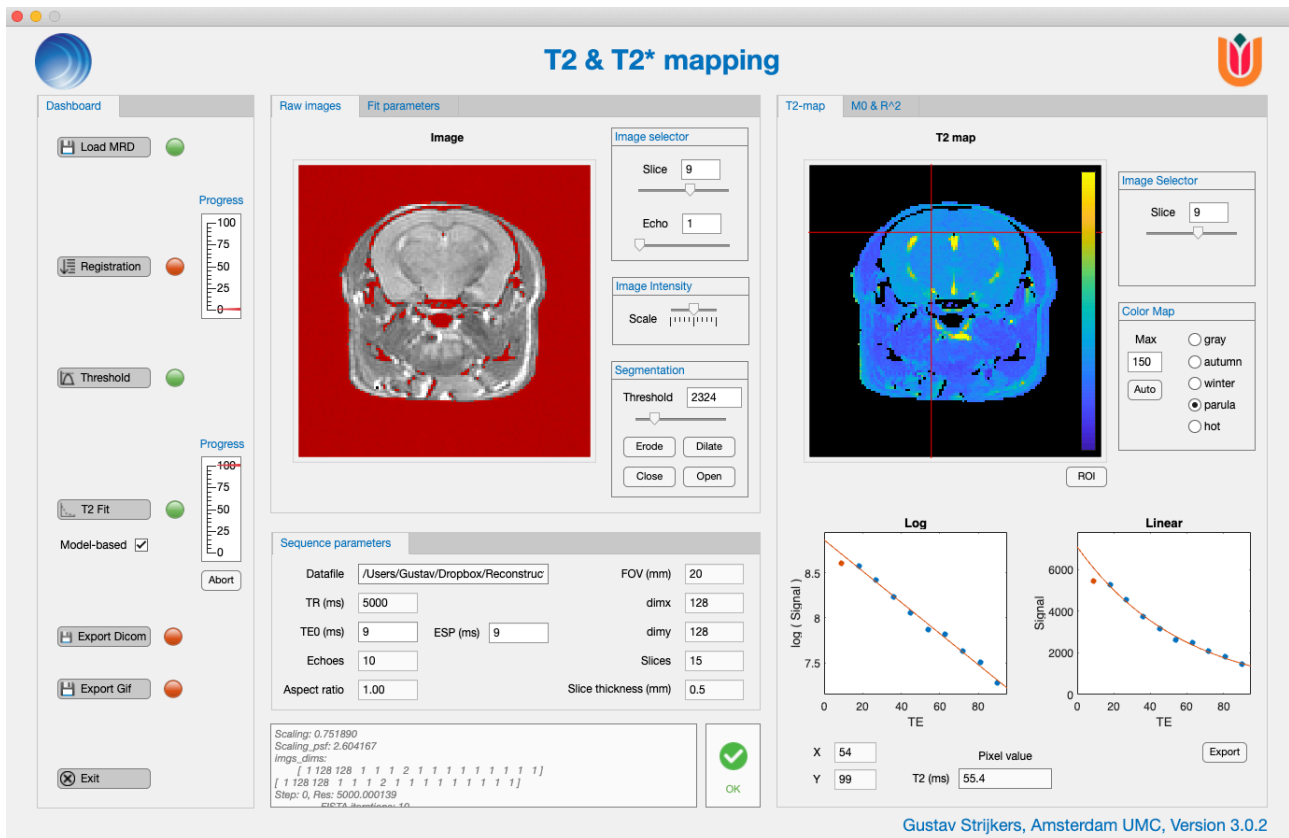

Gustav Strijkers, Amsterdam UMC, Version 3.0.2

### Panel 1: Dashboard

This panel contains the task buttons and parameters that control the reconstruction process. The dashboard tasks need to be completed from top to bottom. A green light next to the task indicates that the task has been completed. Red indicates not completed yet.

### Panel 2: Raw images / Fit parameters

This panel has 2 tabs, raw images and fit parameters.

### Panel 3: Scan parameters

This panel displays the relevant acquisition parameters.

### Panel 4: T2-map / M0 & R<sup>2</sup>

During and after fitting is completed this panel shows the fitted T2-maps. M0-maps and R<sup>2</sup> goodness of fit maps are displayed in the second tab.

### Panel 5: Messages

Displays program status and messages.

## Step 1: Loading data

Press  to import load the MRD raw data file.

T2mapp in needs multi-echo spin-echo data or multi-gradient-echo (FLASH) data.

For multi receiver coil data, select 1 of the MRD raw data files.

Relevant acquisition parameters will be shown in panel 3.

Sequence parameters

Datafile	/Users/Gustav/Dropbox/Reconstruc		FOV (mm)	20
TR (ms)	5000		dimx	128
TE0 (ms)	9	ESP (ms) 9	dimy	128
Echoes	10		Slices	15
Aspect ratio	1.00		Slice thickness (mm)	0.5

## Step 2: Image registration

**NOTE:** this step usually can be omitted.

Use the sliders and/or edit-fields to inspect the slices and images as function of echo-time (TE). The image scale can be adjusted with the Scale slider.

Image

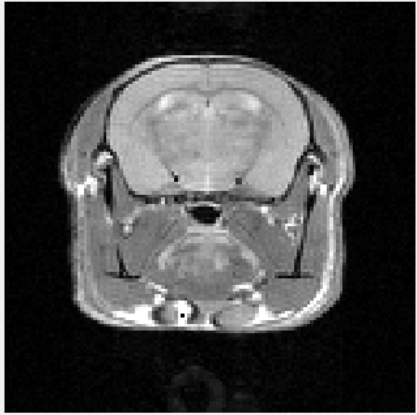


Image selector

Slice

Echo

Image Intensity


Scale

Segmentation

Threshold

Erode Dilate


Close Open

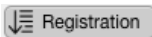In case the images as function of TE are not perfectly registered, press  to perform an affine registration which includes translation, rotation, and scale.




### Step 3: Segmenting the images

To prevent unnecessary fitting of the background pixels, the images need to be segmented first.

Press  Threshold to perform an automatic thresholding segmentation.


The result will look something like the figure below in which the red pixels indicate the regions that will be skipped during T2 fitting.

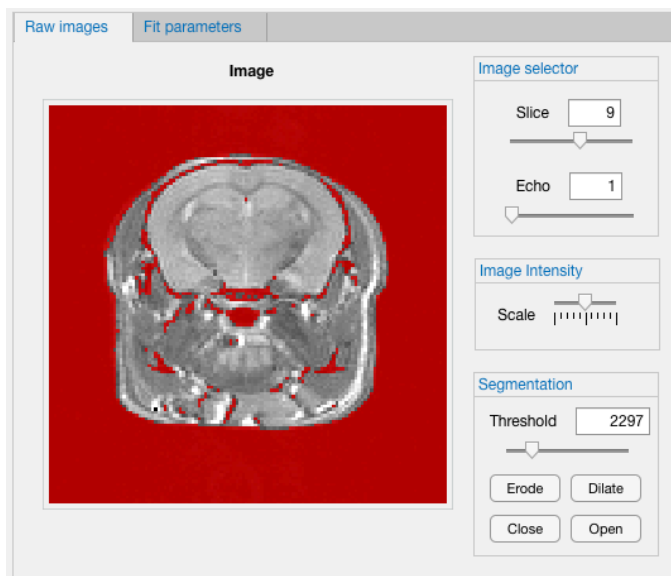

In case automatic thresholding is not satisfying, the threshold can be manually adjusted per slice. The segmentation can be optimized by erode, dilate, close, and open morphological operations.

## Step 4: Fitting T2 maps

Press the Fit parameters tabs to inspect some fit settings.

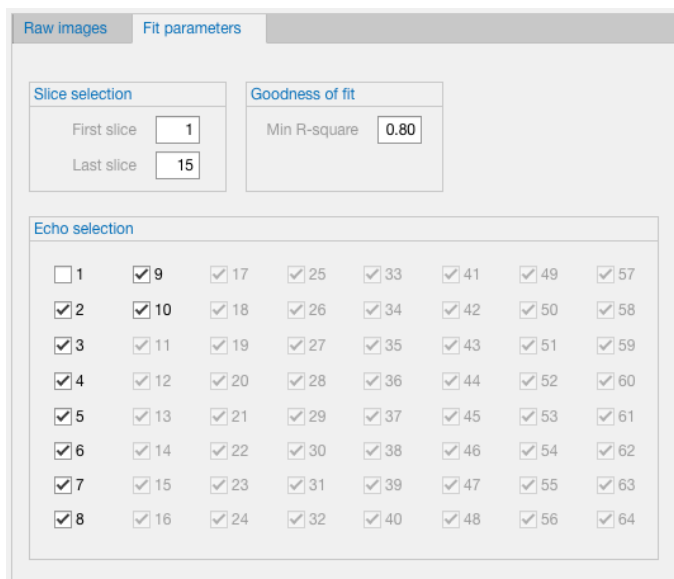

Slice selection		Goodness of fit	
First slice	1	Min R-square	0.80
Last slice	15		

Echo selection							
<input type="checkbox"/> 1	<input checked="" type="checkbox"/> 9	<input checked="" type="checkbox"/> 17	<input checked="" type="checkbox"/> 25	<input checked="" type="checkbox"/> 33	<input checked="" type="checkbox"/> 41	<input checked="" type="checkbox"/> 49	<input checked="" type="checkbox"/> 57
<input checked="" type="checkbox"/> 2	<input checked="" type="checkbox"/> 10	<input checked="" type="checkbox"/> 18	<input checked="" type="checkbox"/> 26	<input checked="" type="checkbox"/> 34	<input checked="" type="checkbox"/> 42	<input checked="" type="checkbox"/> 50	<input checked="" type="checkbox"/> 58
<input checked="" type="checkbox"/> 3	<input checked="" type="checkbox"/> 11	<input checked="" type="checkbox"/> 19	<input checked="" type="checkbox"/> 27	<input checked="" type="checkbox"/> 35	<input checked="" type="checkbox"/> 43	<input checked="" type="checkbox"/> 51	<input checked="" type="checkbox"/> 59
<input checked="" type="checkbox"/> 4	<input checked="" type="checkbox"/> 12	<input checked="" type="checkbox"/> 20	<input checked="" type="checkbox"/> 28	<input checked="" type="checkbox"/> 36	<input checked="" type="checkbox"/> 44	<input checked="" type="checkbox"/> 52	<input checked="" type="checkbox"/> 60
<input checked="" type="checkbox"/> 5	<input checked="" type="checkbox"/> 13	<input checked="" type="checkbox"/> 21	<input checked="" type="checkbox"/> 29	<input checked="" type="checkbox"/> 37	<input checked="" type="checkbox"/> 45	<input checked="" type="checkbox"/> 53	<input checked="" type="checkbox"/> 61
<input checked="" type="checkbox"/> 6	<input checked="" type="checkbox"/> 14	<input checked="" type="checkbox"/> 22	<input checked="" type="checkbox"/> 30	<input checked="" type="checkbox"/> 38	<input checked="" type="checkbox"/> 46	<input checked="" type="checkbox"/> 54	<input checked="" type="checkbox"/> 62
<input checked="" type="checkbox"/> 7	<input checked="" type="checkbox"/> 15	<input checked="" type="checkbox"/> 23	<input checked="" type="checkbox"/> 31	<input checked="" type="checkbox"/> 39	<input checked="" type="checkbox"/> 47	<input checked="" type="checkbox"/> 55	<input checked="" type="checkbox"/> 63
<input checked="" type="checkbox"/> 8	<input checked="" type="checkbox"/> 16	<input checked="" type="checkbox"/> 24	<input checked="" type="checkbox"/> 32	<input checked="" type="checkbox"/> 40	<input checked="" type="checkbox"/> 48	<input checked="" type="checkbox"/> 56	<input checked="" type="checkbox"/> 64

The first and last slice that will be fitted can be selected.

Min R-square is the minimally accepted  $R^2$  value that is accepted after fitting. If the result of a T2 fit in a pixel results in an  $R^2$  value less than Min R-square, then the T2 and M0 values in this pixel will be omitted and set to zero.


The bottom part of this tab shows a list of all the echo numbers. The echos that are selected will be included in the fitting process. This option can be used to exclude certain TEs, for example the first TE, the last TEs below noise level, or in some cases you may only want to fit the even or odd echo numbers.

Press  to start the T2 fitting.

Fitting will be done using a least-squares algorithm according to the equation

$$\text{Signal} = M_0 \exp(-TE/T_2)$$



## Model-based fitting

Tick the model-based checkbox  to perform a model-based reconstruction of the data using the Bart toolbox.

Model-based reconstruction is currently only implemented for T2, not for T2\*.

## Step 5: Fit results

The T2 maps will be shown in the T2-map tab. An M0 map and R<sup>2</sup> map are shown in the second tab. A different color scheme can be chosen for the different maps.

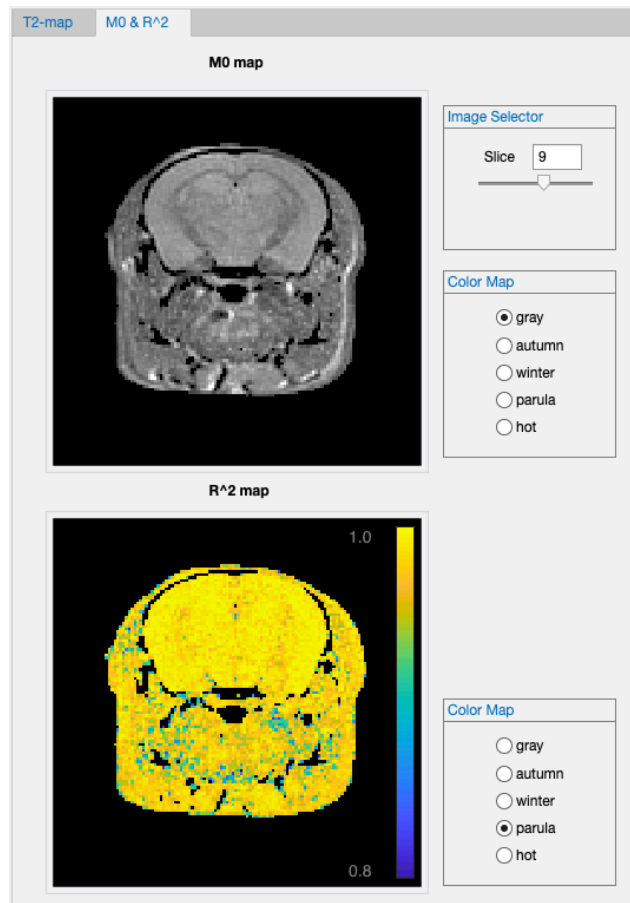
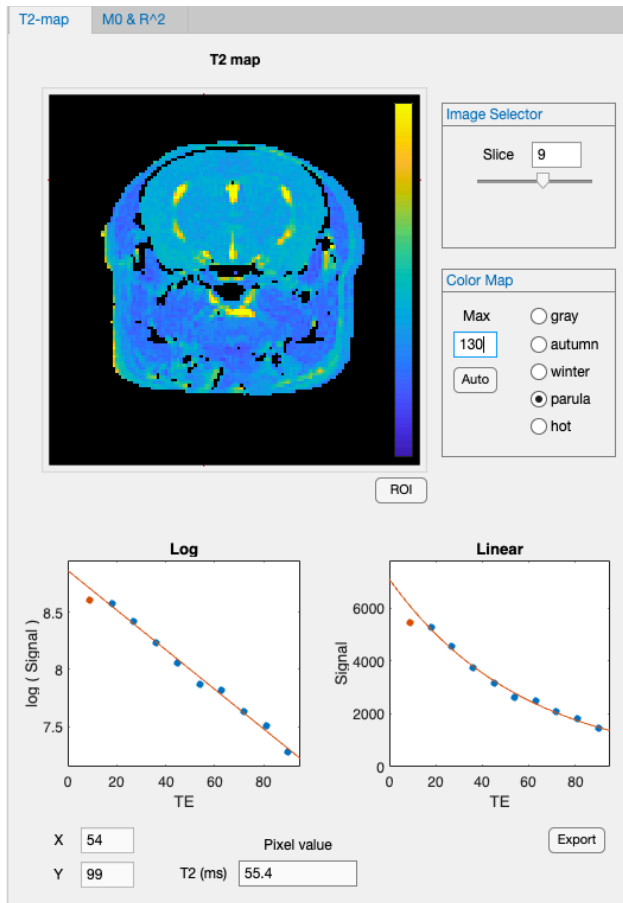




## Step 6: Extracting T2 and T2\* values from pixels and ROIs

When you click on the T2 map with the mouse cursor a cross-hair will be shown. The  $\log(\text{Signal})$  and Signal as function of TE (symbols) and the fit result (red lines) of the cross-hair center pixel, as well as the fitted T2(\*) value of that pixel will be shown.

A red symbol indicates that this TE was omitted in the fitting process, according to the list in the Fit parameters tab (see step 4).

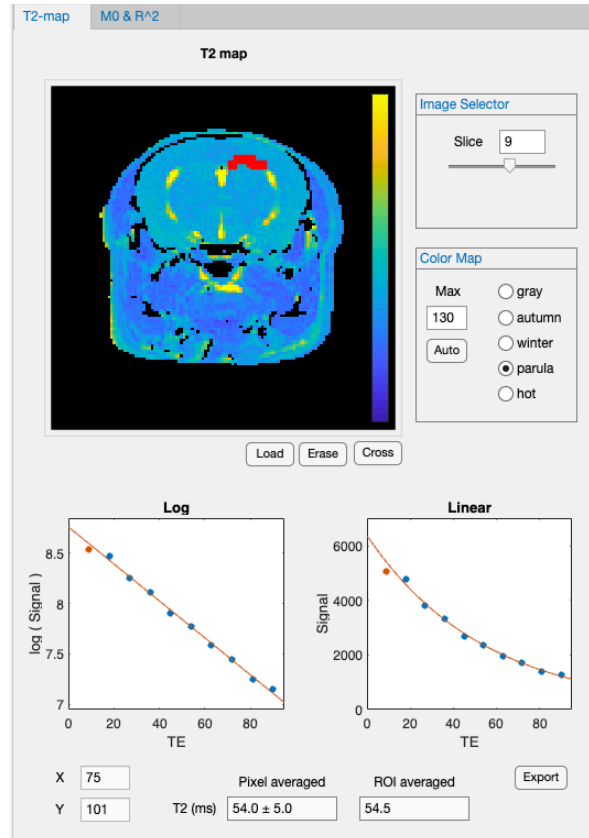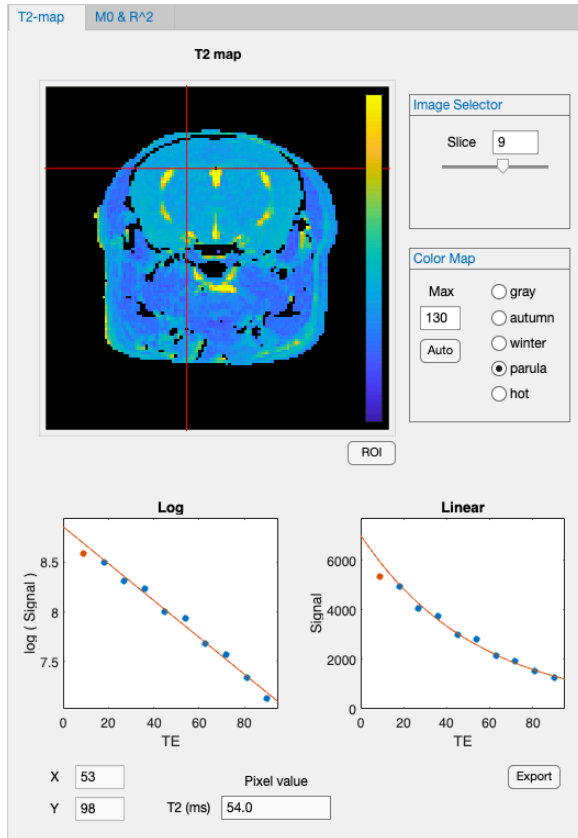

You can toggle between the cross-hair and region-of-interest (ROI) values by pressing the **ROI** button. A single ROI can be 'painted' on the T2 map in one or more slices by mouse clicking on the image.

For the ROI, two T2(\*) values are calculated:

- Pixel averaged: The average of the T2(\*) values in the pixels.  

$$T2_{ROI} = 1 / \text{Sum}(R2_{\text{pixel}})$$
- ROI averaged: The T2(\*) equation is fitted to the average signal intensities in the ROI, producing a single T2(\*) value for that ROI.

The ROI can be removed with **Erase** .

Alternatively you can load a NIFTI ROI file by pressing **Load** .The NIFTI ROI file should have the same matrix dimensions as the T2 map. Such a ROI can be made for example with ITKsnap (<http://www.itksnap.org>) using the corresponding DICOM images.



## Step 7: Exporting the T2 maps

There are two ways to export the T2, M0, and R<sup>2</sup> maps.


### (1) **Export Dicom**

Exports the data in Dicom format for further processing in 3rd party software. The program searches for the Dicom information. If this information is not found, tags will be generated by the program itself. In the latter case the correct image position and orientation information are lost.

### (2) **Export Gif**

Exports the data in Gif format.

## Step 8: Exit

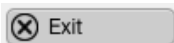Press  to shut down the program.